



Contribution to accurate Spherical Gold Nanoparticles analysis (size, size distribution) by SpICPMS and SAXS

V Geersten, O. Taché, E Barruet, J L Lacour, F. Gobeaux, Antoine Thill

► To cite this version:

V Geersten, O. Taché, E Barruet, J L Lacour, F. Gobeaux, et al.. Contribution to accurate Spherical Gold Nanoparticles analysis (size, size distribution) by SpICPMS and SAXS. Workshop on Reference Nanomaterials: Current situation, measurement methods, standardization and future needs, May 2018, Berlin, Germany. 10.1107/S1600576716012127 . cea-02328404

HAL Id: cea-02328404

<https://cea.hal.science/cea-02328404>

Submitted on 23 Oct 2019

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

Contribution to accurate Spherical Gold Nanoparticles analysis (size, size distribution) by spICPMS and SAXS


Dr. V. Geersten, O. Taché, E. Barruet, J.L. Lacour, Dr. F. Gobeaux and Dr. A. Thill

NIMBE (Nanoscience & Innovation Laboratory) / Fundamental Research Division / CEA Paris-Saclay
91191 Gif sur Yvette Cedex FRANCE

Small-Angles X-Ray Scattering (SAXS) has been established as a metrological method for the determination of nanoparticles size and size distribution. Modern SAXS Laboratory experiments¹, by involving synchrotron-based instrumentation at lower price and very stable X-ray source, are more and more used in nanomaterials domain. In the frame of the EMPIR Innanopart project, we have developed a methodology for the size, size distribution and concentration determination of spherical nanoparticles. This protocol involves a precise sample preparation, and a set of homemade software tools for the data processing - from the acquisition, the absolute scaling, to the analysis.

spICPMs is not a metrological traceable technique but has many strengths to become a useful complement of nanoparticle characterization methods such as SAXS and microscopy. It can also measure highly diluted nanoparticles suspensions which is not the case of Dynamic Light scattering (DLS) or SAXS. Finally, ICPMS analyzes inorganic ions in liquid solution in a very large range of concentration, which should allow linear diameter measurement range over at least 3 orders of magnitude.

In this work, we confront spICPMS with SAXS in order to investigate the method and the developed protocols on a set of commercial spherical Gold Nanoparticles.

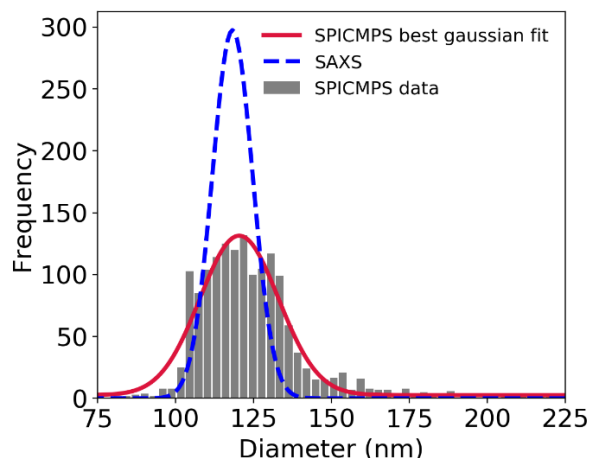

Comparison between SAXS and spICPMS method
for the determination of size of spherical Gold
Nanoparticles

¹ Olivier Taché et al., « MOMAC: a SAXS/WAXS laboratory instrument dedicated to nanomaterials », *Journal of Applied Crystallography* 49, n° 5 (1 octobre 2016): 1624- 31, <https://doi.org/10.1107/S1600576716012127>.