

# Inductively-Coupled Microcoils and Solution Flow for Increased NMR Sensitivity

Guillaume Carret, Thomas Berthelot, Céline Boutin, Patrick Berthault

### ▶ To cite this version:

Guillaume Carret, Thomas Berthelot, Céline Boutin, Patrick Berthault. Inductively-Coupled Microcoils and Solution Flow for Increased NMR Sensitivity. 60th Experimental Nuclear Magnetic Resonance Conference, Apr 2019, Pacific Grove, United States. cea-02329455

## HAL Id: cea-02329455 https://cea.hal.science/cea-02329455v1

Submitted on 23 Oct 2019

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers. L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

#### Inductively-Coupled Microcoils and Solution Flow for Increased NMR Sensitivity


#### Guillaume Carret, Thomas Berthelot, Céline Boutin, Patrick Berthault

#### NIMBE, CEA, CNRS, Paris-Saclay University, CEA Saclay, 91191 Gif-sur-Yvette, France

To circumvent the sensitivity problems inherent to slow return of magnetization to equilibrium in liquid-state NMR we present a 3D-printed NMR device based on a mini bubble-pump associated with fluidics and micro-detection that can be installed in every commercial liquid-state NMR probehead. It is based on the use of a closed-loop circuit of the solution near the NMR magnetic center, which presents two main advantages: i) pre-polarization is achieved for the whole solution volume, ii) this volume can be reduced to tens of microliters.

This device (dubbed WIFI-NMRS) is inserted from the top of the NMR magnet, as would be the case for any liquid sample. A gas flow driven by a programmable syringe pump actuates a mini bubble-pump which leads to circulation of

Vernier.

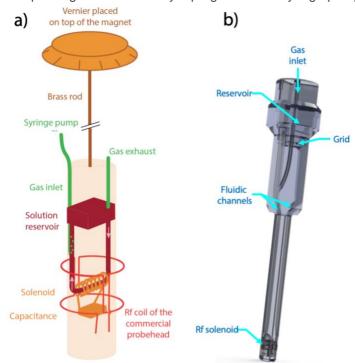

*Figure 1.* a) Principle of the WIFI-NMRS device b) 3D rendering of the NMR insert (8mm-version).

The inductive devices also benefit from all the capabilities of the host probe, allowing heteronuclear experiments even with a mono-tuned microcoil: <sup>1</sup>H decoupling (Figure 2), X-<sup>1</sup>H experiments. Furthermore, the coupling of the microcoil to one of the host probe channels creates two resonant frequencies which can be modulated to extend the probe capabilities, either to reach frequencies unattainable with the host probe alone, or to use both of these resonance frequencies to study two different nuclei with the optimized detection allowed by the microcoil. To increase the signal-to-noise ratio one step further, this system can also be used to efficiently dispense gaseous species such as hyperpolarized xenon and parahydrogen to the solution.

static on flow 200 150 100 50 ppm

the liquid sample. A part of this circuit crosses the NMR

detection region consisting of a micro-solenoid. A brass rod fixed on the upper part of the insert and ended by a

Vernier placed on top of the magnet enables angular

positioning of the micro-coil (Figure 1). The NMR insert has the shape of an NMR tube with its spinner, in order to

be easily positioned inside the probehead. Only the diameter of the lower part varies to be adjusted to the

commercial probehead (from 5 to 20 mm). Optimization of

the coupling between the two coils is made by manually

adjusting the angular position of the insert thanks to the

As an example of the performance of the approach for


slowly relaxing nuclei, magnetic resonance spectroscopy

and imaging experiments have been performed with various samples and nuclei (<sup>13</sup>C, but also <sup>7</sup>Li, <sup>89</sup>Y etc.). The approach of using a circular flow has also a great interest for the continuous monitoring of chemical reactions, and particularly those involving modification of

guaternary carbons. The flow motion homogenizes the

reactant mixture, and as detection of poorly relaxing

nuclei is facilitated, more quantitative results are obtained.

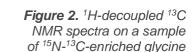

Finally since the quality factor depends on the host probe, the device will be very even more powerful when used in conjunction with cryoprobes.

Reference. G. Carret, T. Berthelot, P. Berthault, Analytical Chemistry 90 (2018) 11169–11173, DOI: 10.1021/acs.analchem.8b01775

Acknowledgement. Support from the French Ministry of Research (ANR 17-LCV2-0002-01 LabCom DESIR) is acknowledged.