

Qualitative and quantitative analysis of materials by LIBS in nuclear environment how multivariate methods can help.

J-B. Sirven, M. El Rakwe, N. Coulon, D. L'Hermite, E. Vors

► To cite this version:

J-B. Sirven, M. El Rakwe, N. Coulon, D. L'Hermite, E. Vors. Qualitative and quantitative analysis of materials by LIBS in nuclear environment how multivariate methods can help.. LIBS 2016, Sep 2015, Chamonix-Mont-Blanc, France. cea-02439451

HAL Id: cea-02439451 https://cea.hal.science/cea-02439451v1

Submitted on 26 Feb 2020 $\,$

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers. L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Ceaden





QUALITATIVE AND QUANTITATIVE ANALYSIS OF MATERIALS BY LIBS IN NUCLEAR ENVIRONMENT: HOW MULTIVARIATE METHODS CAN HELP

J.-B. Sirven, M. El Rakwe, N. Coulon, D. L'Hermite, E. Vors

CEA, Nuclear Energy Division, Department of Physical Chemistry F-91191 Gif sur Yvette, France

jean-baptiste.sirven@cea.fr

LIBS 2016, Sept. 12 – 16, 2015, Chamonix-Mont-Blanc

Ceaden NUCLEAR ENVIRONMENT





FROM RESEARCH TO INDUSTRY

Ceaden NUCLEAR ENVIRONMENT





FROM RESEARCH TO INDUSTRY

Ceaden NUCLEAR ENVIRONMENT








Ceaden NUCLEAR ENVIRONMENT










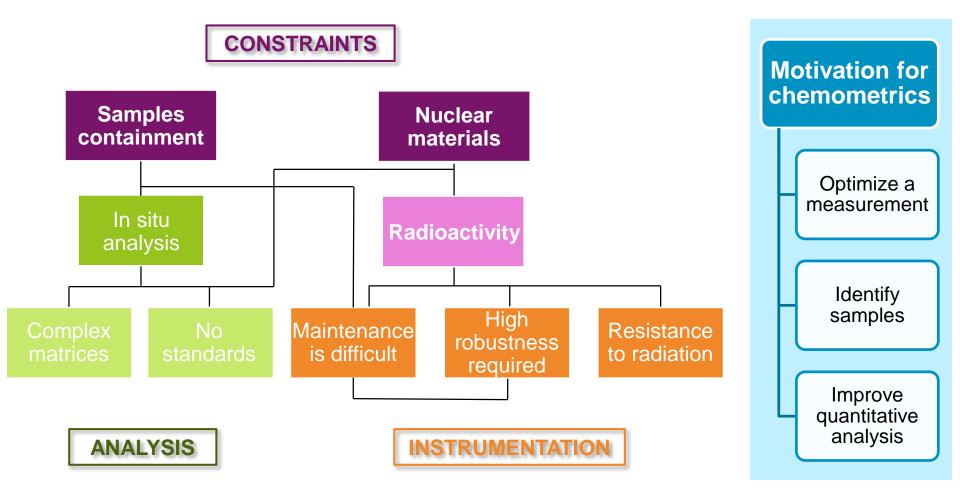

1/3

EXPERIMENTAL DESIGNS TO IMPROVE THE ROBUSTNESS OF THE ANALYSIS

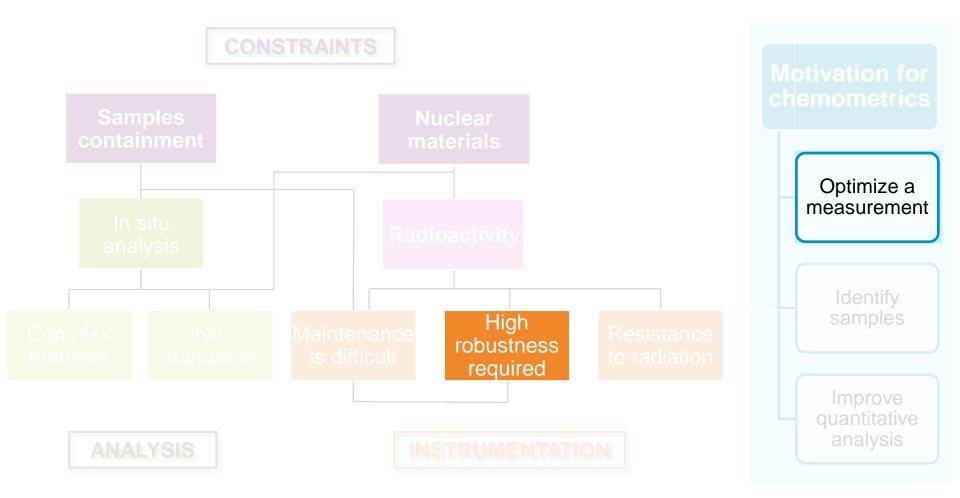

is non-linear δ td Factors are strongly **coupled** sample

Objective: To optimize experimental parameters in ordre to obtain the best analytical performances, with a user point of view and a minimum of experimental trials

OPTIMIZATION OF INSTRUMENTAL


Ceaden PARAMETERS USING A DESIGN OF

detection

gate delay (t_d)

Φ

EXPERIMENTS

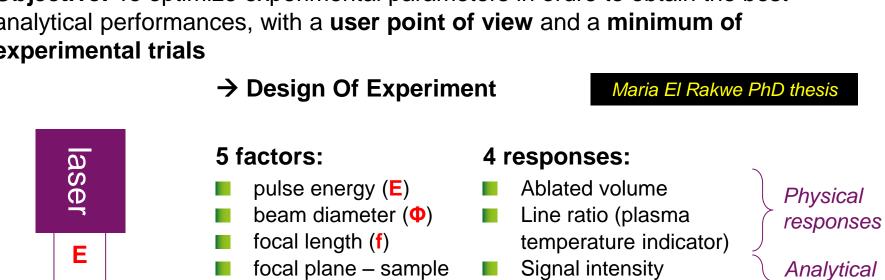

Taking into account that system The relation between factors and responses

Signal repeatability

responses

Sirven et al. | Qualitative and quantitative analysis of materials by LIBS in nuclear environment: how multivariate methods can help | PAGE 9

Experiment:

266 nm Nd:YAG laser Czerny-Turner monochromator, 300 mm, 2400 g/mm grating ICCD


Sample:


SS59 standard: Fe 99.6 % ; Mn 0.12 %

Central composite DOE (32 experiments)

- 1. Modeling of relation between factors and responses
- 2. Simultaneous optimization of analytical responses
- 3. Check of experimental reproducibility
- 4. Experimental validation of the predictive ability of the model



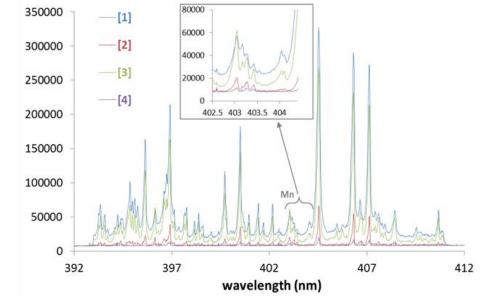

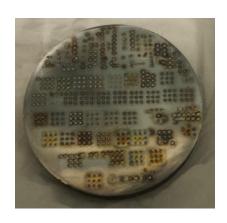






OPTIMIZATION OF INSTRUMENTALCeaderPARAMETERS USING A DESIGN OFEXPERIMENTS

Experimental reproducibility

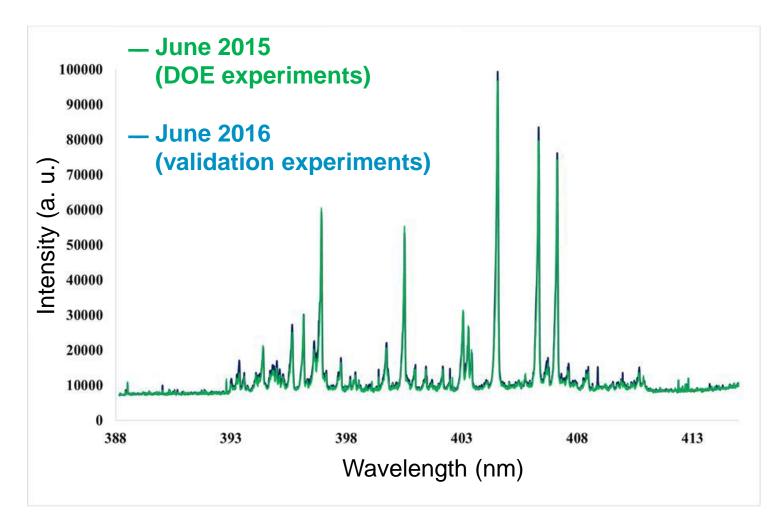





Model validation


Responses	Predicted values	Experimental values	Relative deviation
Mn signal intensity	7260 ± 1740	7530	- 3.7%
Mn signal RSD	1 % (< 6 %)	3.5 %	+ 250%
Ablated volume (µm ³)	(1.45 ± 0.16) 10 ⁵	1.27 10 ⁵	- 12%
Line ratio (plasma temperature indicator)	1.7 ± 1.7	2.9	+ 68%

The model is satisfactorily predictive → it can be used:

- To determine the effect of the chosen experimental parameters on different responses (e.g. SBR, LOD...)
- To estimate the **instrumental robustness** in a given configuration



MATERIALS IDENTIFICATION BY ONSITE LIBS SYSTEMS

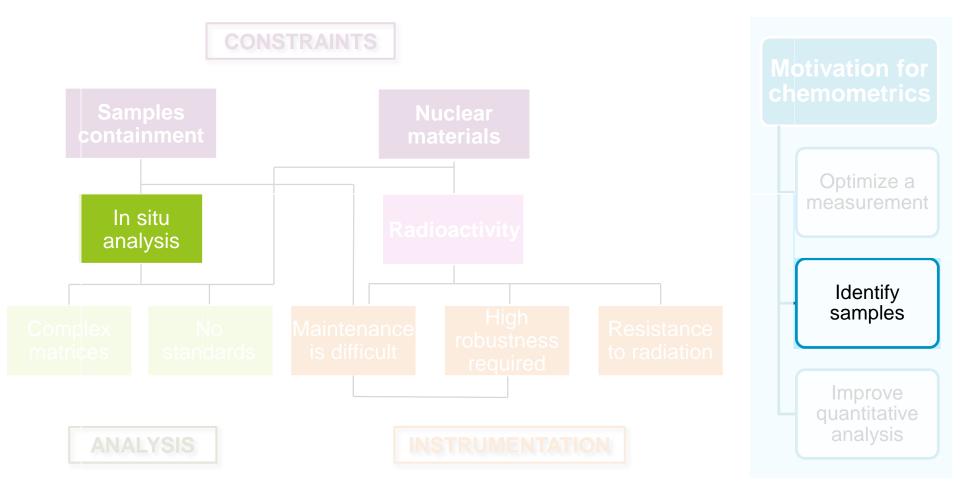

Inventory before nuclear decommissioning Control of contamination Control of contamination

Identification of wastes prior to treatment / sorting

2 technologies:

Portable

Fibered

Motivation:






RICA robot developed at CEA for in situ measurements in nuclear environment





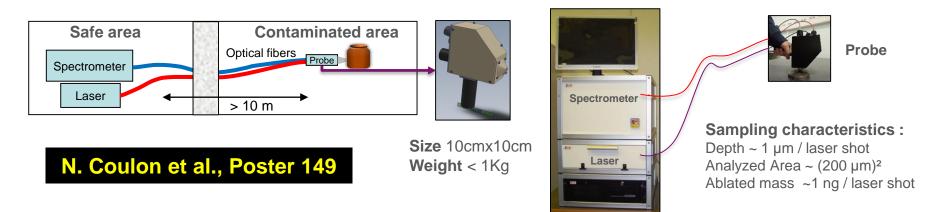





Ceaden ONSITE IDENTIFICATION OF MATERIALS


15500-15000-14500-14000-13500-13000-12500-12000-11500-


11000




Examples of data treatment

Single-shot, low resolution spectra (portable system)

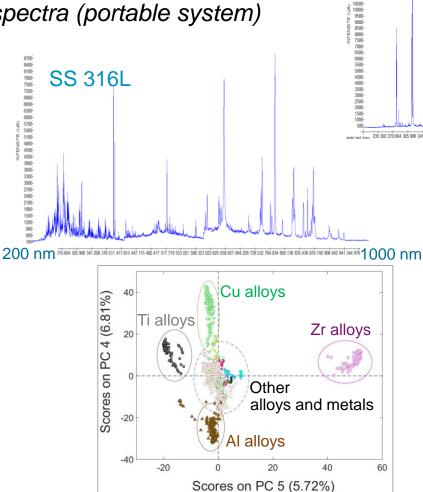


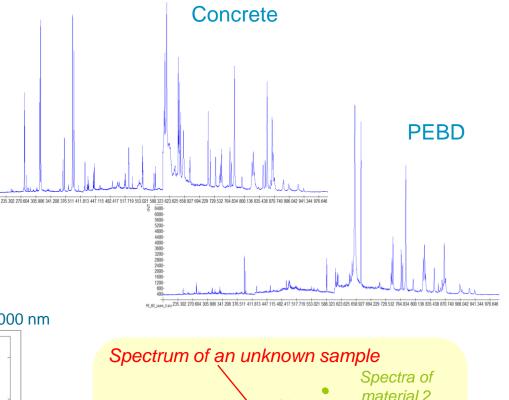


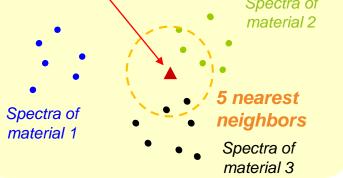

ONSITE IDENTIFICATION OF MATERIALS



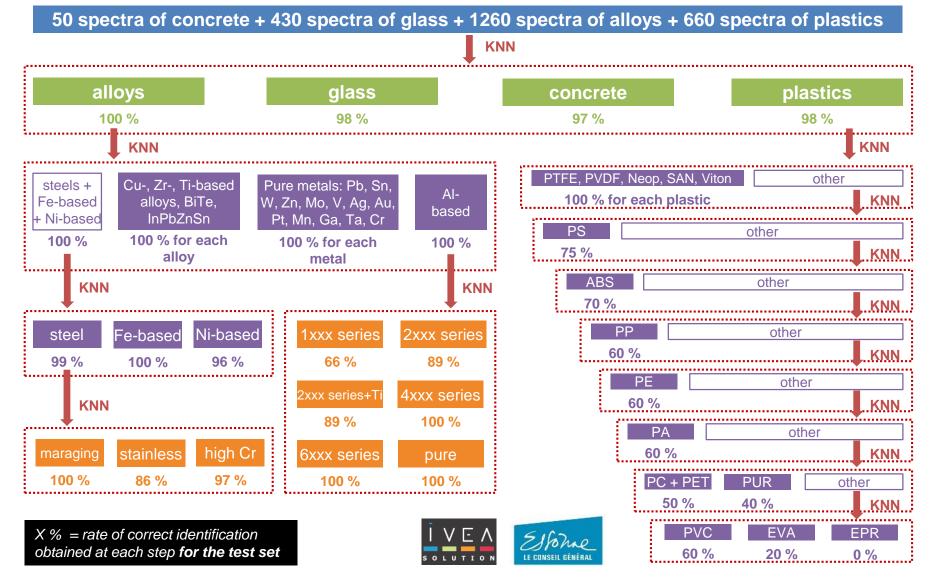

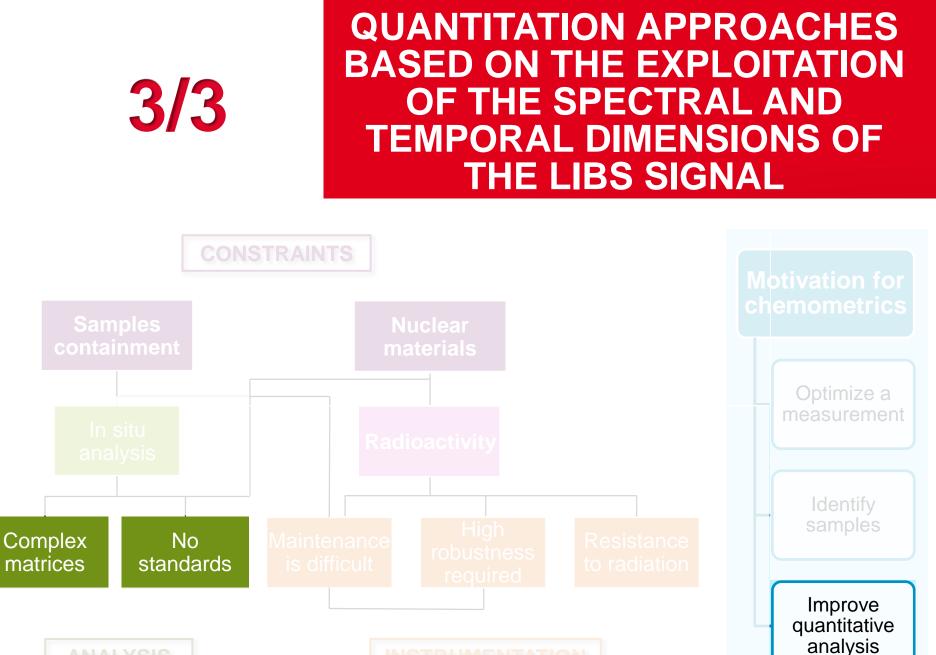

ANALYSIS

INSTRUMENTATION

FROM RESEARCH TO INDUSTRY

Ceaden calibration by LIBS




Calibration of [Ni] in...



• Zr matrix


Cu matrix

- Gate delay = $1.5 \ \mu s$
- Gate width = $3 \mu s$

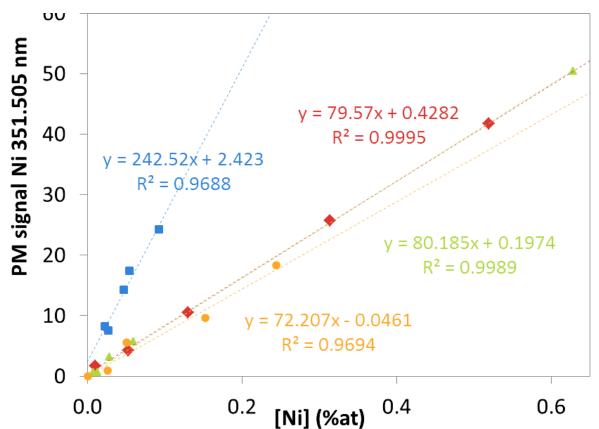

Ceaden EFFECT OF TIME RESOLUTION ON CALIBRATION LINES



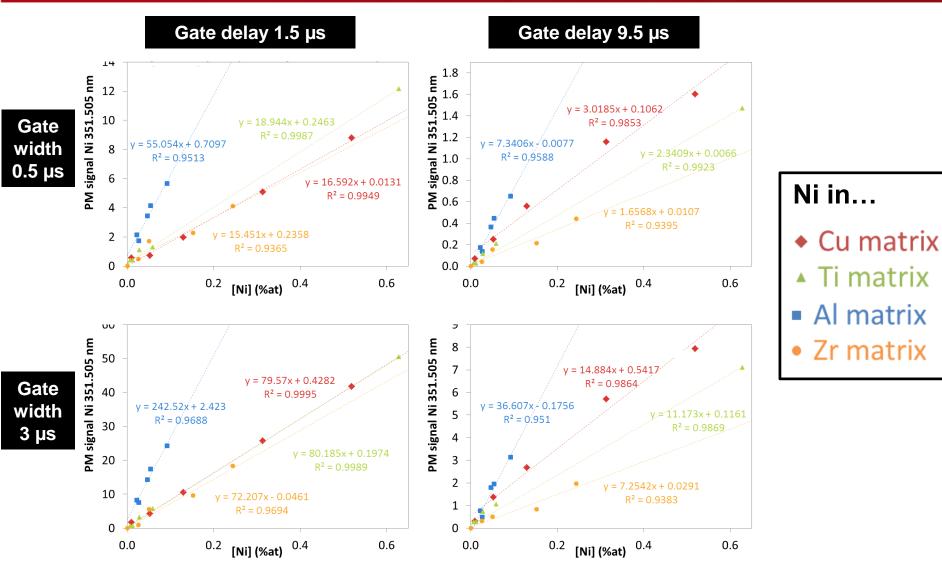


Sirven et al. | Qualitative and quantitative analysis of materials by LIBS in nuclear environment: how multivariate methods can help | PAGE 18

Ceaden CALIBRATION USING TIME-RESOLVED SPECTRA




What if we take into account the temporal dimension of the LIBS signal?

- 1. Does it improve analytical performances?
- 2. Can we take advantage of it to correct matrix effects?

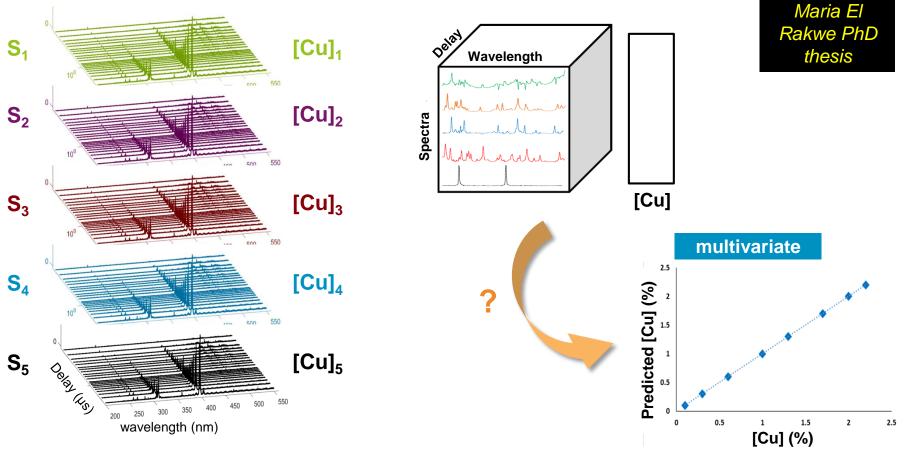


Ceaden A PRACTICAL EXAMPLE



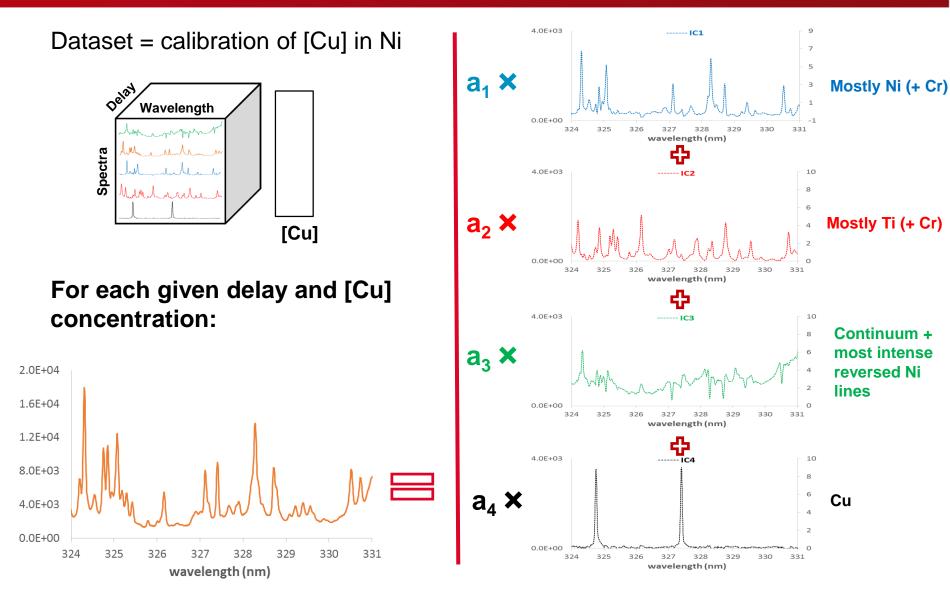

CALIBRATION LINE VS PLS VS ICA-PLS




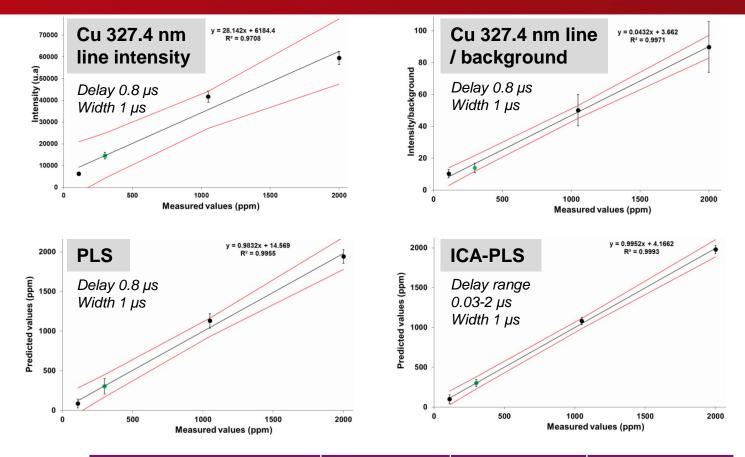

	Calibration	Trueness	Precision	LOD (ppm)	
	Cu line intensity			900	
	Cu line / background			400	
	PLS			350	
Sirven et al	ICA-PLS			180	can help PAGE 21

Ceaden MULTI-MATERIALS QUANTITATIVE ANALYSIS



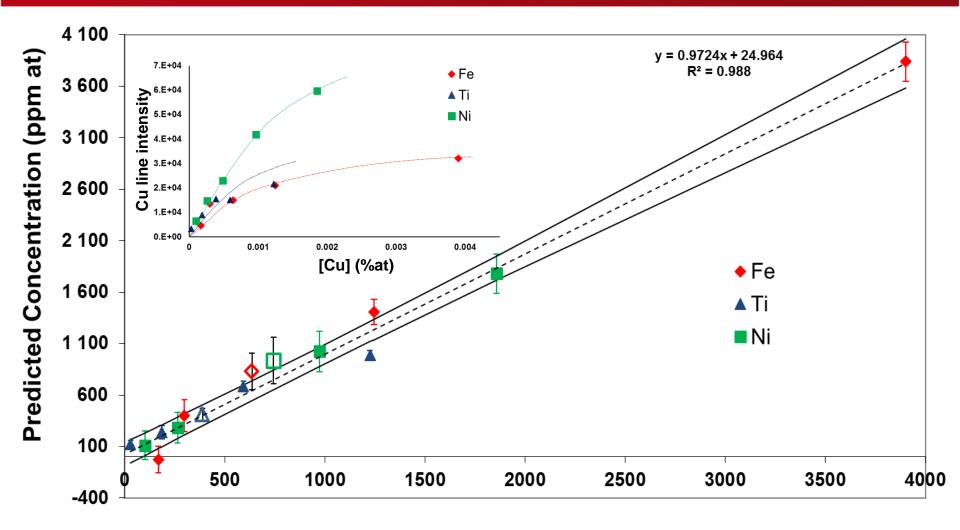

Measured Concentration (ppm at)


Ceaden conclusion



Elemental analysis in nuclear environment:

- Strong instrumental constraints
- Complex spectra
- Few (sometimes no) standards available
- Approaches based on Boltzmann modeling (CF-LIBS, Cσ, FP...) difficult to implement

Chemometrics is useful to overcome (to a certain extent) those difficulties (but it is not magic... and requires careful validation...)

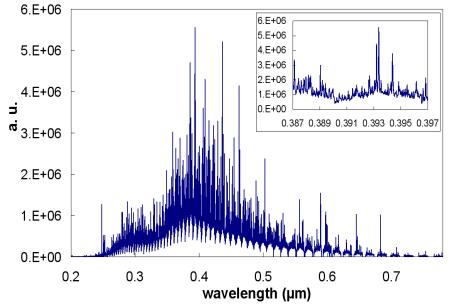

Sirven et al. | Qualitative and quantitative analysis of materials by LIBS in nuclear environment: how multivariate methods can help | PAGE 23

Acknowledgements:

Ayoub Aouam Rémi Beaubaton Alix Dehayem-Massop Guillaume Gallou Tiphanie Lelong Gilles Moutiers Laurie Pontreau Sébastien Robino Douglas N. Rutledge Thomas Vercouter

Thank you for your attention!

jean-baptiste.sirven@cea.fr

Commissariat à l'énergie atomique et aux énergies alternatives	DEN
Centre de Saclay 91191 Gif-sur-Yvette Cedex	DANS
T. +33 (0)1 69 08 43 71	DPC
	SEARS
Etablissement public à caractère industriel et commercial RCS Paris B 775 685 019	LANIE

eaden Abstract



Elemental analysis in the nuclear field is essential in a number of situations, from the frontend of the fuel cycle to the waste management, either for process control, compliance checking, or waste characterization prior to conditioning. Yet, it is special for two main reasons.

First, due to their radioactivity, the samples are confined inside containments such as controlled areas, glove boxes or hot cells. Those safety and security constraints inherent in nuclear environment strongly restrict the possibility of sampling, transportation and handling of materials. Therefore, onsite or even in situ analytical techniques, such as LIBS, are of high interest, as they enable to limit the workers exposure to radiation and contamination, as well as the production of waste and effluents generated by conventional laboratory analysis. One important consequence is that the instruments must be very robust and require as little maintenance as possible.

Secondly, the matrices can be very complex due to the presence of many interfering elements including actinides and lanthanides. And a large variety of materials need to be characterized (fuel, alloys, glasses, concrete, solutions, etc.), both qualitatively and quantitatively. In this latter case, calibration samples are not always available and alternative approaches must be developed.

As shown by the abundant literature on the subject in other areas of analytical spectroscopy, multivariate methods are very helpful to cope with such issues. We will show how similar strategies can be implemented in LIBS, taking into account the peculiarities of the technique. Different applications of interest in the nuclear field will be reviewed: materials identification by onsite LIBS systems; use of experimental designs to improve the robustness of the analysis; innovative quantitation approaches based on the exploitation of the spectral and temporal dimensions of the LIBS signal; and isotopic analysis by an original method based on the self-absorption of intense emission lines.