



**HAL**  
open science

## Development and comparison of two high accuracy methods for uranium concentration in nuclear fuel: ID-TIMS and K-edge densitometry

Alexandre Quemet, Alexandre Ruas, Eric Esbelin, Vincent Dalier, Cédric Rivier

### ► To cite this version:

Alexandre Quemet, Alexandre Ruas, Eric Esbelin, Vincent Dalier, Cédric Rivier. Development and comparison of two high accuracy methods for uranium concentration in nuclear fuel: ID-TIMS and K-edge densitometry. *Journal of Radioanalytical and Nuclear Chemistry*, 2019, 321 (3), pp.997-1004. 10.1007/s10967-019-06670-y . cea-03470768

**HAL Id: cea-03470768**

**<https://cea.hal.science/cea-03470768v1>**

Submitted on 8 Dec 2021

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

1

## **Title**

2 Names of the authors: Alexandre Quemet<sup>1</sup>, Alexandre Ruas<sup>2</sup>, Eric Esbelin<sup>1</sup>, Vincent  
3 Dalier<sup>1</sup> and Cédric Rivier<sup>1</sup>

4 Title: Development and comparison of two high accuracy methods for uranium  
5 concentration in nuclear fuel: ID-TIMS and K-edge densitometry

6 Affiliation(s) and address(es) of the author(s):

7 - <sup>1</sup> CEA, DEN, DMRC, SA2I, L2AT, BP17171, 30207 Bagnols sur Cèze, France

8 - <sup>2</sup> Onsite Laboratory Team, Nuclear Material Laboratory, Office of Safeguards  
9 Analytical Services, Department of Safeguards, International Atomic Energy  
10 Agency, Tokyo Regional Office, Seibunkan Bldg, 9F, 1-5-9 Iidabashi, Chiyoda-ku,  
11 Tokyo 102-0072, Japan

12 E-mail address of the corresponding author: alexandre.quemet@cea.fr

13



36        **Introduction**

37        Uranium detailed accountancy is of prime interest in the nuclear industry, at various stages  
38        of the nuclear fuel cycle and for safeguards purposes [1–3]. ATALANTE is a nuclear  
39        facility of the French Alternative Energies and Atomic Energy Commission (CEA) and  
40        dedicated to the research on the nuclear spent fuel reprocessing. Studies such as nuclear  
41        fuel dissolution, actinides separation and precipitation are performed, requiring to  
42        determine the uranium concentration with high accuracy.

43        Among the specific analytical techniques available in such facility, K-edge densitometry  
44        (KED) is an effective technique to determine uranium assay in nuclear fuel dissolution  
45        solution [4, 5]. KED is an X-ray absorbance spectrometry technique based on the Beer-  
46        Lambert's law. The uranium concentration in solution is proportional to the X-ray  
47        attenuation. High uranium concentrations, typically above 50 g L<sup>-1</sup>, can be measured with  
48        high accuracy. For lower concentration, spectral discontinuities make the spectrum process  
49        difficult and uncertainties increase significantly. K-edge spectrometry is a nondestructive,  
50        non-invasive, relatively fast technique and requires little sample preparation.

51        Isotope Dilution (ID) is a well-known method for quantification based on internal  
52        calibration [6, 7]. The ID principle is to mix a known amount of sample containing a known  
53        isotope composition with a spike solution containing a known concentration and isotope  
54        composition. The spike solution contains the same analyte as the sample with a different  
55        isotope composition. The isotope ratio of the sample – spike mixture reflects the sample  
56        analyte concentration. The ID methodology requires only isotope ratio measurements and  
57        knowing masses involved in the sample – spike mixture. As it is based on isotope ratio  
58        determination, an accurate technique for isotope ratio measurements must be preferred.  
59        Thermal-Ionization Mass Spectrometry (TIMS) is the reference technique for the accurate  
60        determination of uranium isotope ratios [8–13]. The ID methodology coupled to TIMS for  
61        the isotope ratio measurements (ID-TIMS) is then a powerful combination for the uranium  
62        concentration determination.

63 The present study was performed in the context of the “2017 Nuclear Material Round  
64 Robin” organized by the International Atomic Energy Agency (IAEA) which aims at  
65 determining the uranium isotope ratio and mass fraction in nuclear materials. The isotope  
66 ratios determination was presented in a previous paper [10]. This paper presents a  
67 comparison between the ID-TIMS and the K-edge densitometry for the uranium  
68 concentration determination with high accuracy in a uranium pellet.

## 69 **Experimental**

### 70 **Materials, reagents and certified reference materials**

71 All solutions were prepared in polypropylene vials, except for uranium solutions, which  
72 were prepared in PFA vials. 1 mol L<sup>-1</sup>, 3 mol L<sup>-1</sup> and 8 mol L<sup>-1</sup> nitric acid solutions were  
73 prepared by diluting high purity nitric acid (Merck, Suprapur) in deionized water  
74 (resistivity: 18 2 MΩ.cm). A high accurated scale (Mettler-Toledo, WXTF 205) was used  
75 to prepare all solutions.

76 The spike solution used in the ID methodology was the IRMM 054 certified reference  
77 material (CRM) provided by the Joint Research Centre of the European Commission  
78 (EC-JRC). This CRM contains high amounts of the <sup>235</sup>U isotope and is certified for the  
79 concentration and the isotope ratios. This solution was diluted gravimetrically in order to  
80 obtain a uranium concentration near 10 µg g<sup>-1</sup>.

### 81 **Sample preparation**

82 The inter-laboratory comparison (ILC) sample is a low enrichment uranium oxide pellet.  
83 The K-edge densitometer available in the laboratory and the TIMS require working with  
84 solutions. The ILC uranium pellet was dissolved in hot nitric acid (8 mol L<sup>-1</sup>). This solution  
85 (pellet + nitric acid) was heated to 135 °C in a PFA vial until complete dissolution. The  
86 dissolution solution had about 250 g L<sup>-1</sup> uranium concentration [10].

87 Instrumental

88 *Thermal Ionization Mass Spectrometer*

89 The Thermo Fisher Triton and VG Sector 54 TIMS used for the experiments and the  
90 deposit technique were previously described in detail [14]. The measurements were  
91 performed using Faraday cups coupled to  $10^{11} \Omega$  current amplifiers (hereafter referred to  
92 as FC 11). All the isotope ratio measurements were performed with the total evaporation  
93 method (TE method) described previously [10]. The  $^{238}\text{U}^+$  ion beam target intensity was  
94 set at 1 V and 2.2 V for the VG Sector 54 and for the Triton TIMS, respectively.

95 The uranium concentration ( $[\text{U}]_S$ ) was determined using isotope dilution according to Eq.  
96 (1):

$$[\text{U}]_S = [\text{U}]_T \cdot \frac{m_T}{m_S} \cdot \frac{M_S}{M_T} \cdot \frac{(235)_T}{(238)_S} \cdot \frac{R_{\text{Mix}} - R_T}{1 - R_{\text{Mix}} \cdot R_S} \quad (1)$$

97 Where  $T$  refers to the spike (or tracer),  $Mix$  refers to the sample - spike mixture and  $S$  refers  
98 to the sample;  $m$  are the masses involved in the mixture,  $M$  are the atomic weights,  $(235)$   
99 and  $(238)$  are the isotope abundances.  $R_{mix}$  and  $R_T$  correspond to the  $^{238}\text{U}/^{235}\text{U}$  mixture and  
100 spike isotope ratio, respectively.  $R_S$  corresponds to the  $^{235}\text{U}/^{238}\text{U}$  sample isotope ratio. The  
101 isotopic composition of the pellet (needed for the  $[\text{U}]_S$  calculation by ID-TIMS) was  
102 determined with high accuracy measurements in the previous study [10].

103 The uranium mass fraction determination by ID-TIMS was performed using the two TIMS  
104 available in the laboratory with 3 different determinations. First, 3 diluted dissolution  
105 solutions ( $[\text{U}] \approx 10 \mu\text{g g}^{-1}$  and hereafter referred to as diluted solution 1, 2 or 3) were  
106 prepared gravimetrically from the initial dissolution solution ( $[\text{U}]_i \approx 250 \text{ g.L}^{-1}$ ). Then, 5  
107 (diluted solution 1 – spike) mixtures, 3 (diluted solution 2 – spike) mixtures and 2 (diluted  
108 solution 3 – spike) mixtures were prepared gravimetrically. Each mixture was evaporated  
109 after preparation, and then dissolved again in  $10 \mu\text{L}$  of  $1 \text{ mol L}^{-1}$  nitric acid.  $1 \mu\text{L}$  (about  
110  $400 \text{ ng}$  of U) of this final solution was deposited onto a filament for each individual

111 measurement. The  $^{238}\text{U}/^{235}\text{U}$  isotope ratio of each mixture was analyzed using the TE  
112 method. The  $^{238}\text{U}/^{235}\text{U}$  ratios ranged from 1.08 to 1.23. The 5 (diluted solution 1 – spike)  
113 mixtures were analyzed using the VG Sector 54 TIMS. The uranium concentration was the  
114 average of the 5 individual measurements. The (diluted solution 2 and 3 – spike) mixtures  
115 were analyzed using the Triton TIMS and the respective uranium concentrations were  
116 determined by computing the average of each respective individual measurements. Finally,  
117 the uranium mass fraction of the pellet was the average of the 3 groups of determination.

### 118 *K-edge densitometry*

119 The nuclearized KED used for the experiments was previously described in details [4]. The  
120 X-ray generator and X-ray tube (W anode, MG165) are marketed by Yxlon. The tube  
121 configuration used in the laboratory is 150 kV and 10 mA for assay sample measurements,  
122 or 5 mA for blank reference solution measurements. The detectors are hyper-pure Ge  
123 crystals (Canberra EGX 200-10, active area 200 mm<sup>2</sup>) cooled with liquid nitrogen. The  
124 acquisition electronics are Canberra DSA-1000 coupled with Genie 2000 software. The  
125 spectral data were processed using a CEA developed software using standardless  
126 algorithms [4]. Measurements are performed without any calibration except for energy.  
127 The typical counting time was 1000 s.

128 This device is located in the back area of the analysis hot cell line of ATALANTE. The  
129 sample, packaged in small polypropylene vials, is transferred *via* a pneumatic line from the  
130 hot cell to the KED. The minimum volume required to ensure the X-rays are absorbed by  
131 the sample in this system is about 1.5 mL.

132 The dissolution solution was analyzed by KED without any further preparation except vial  
133 transfer. The result was obtained from the average of 5 individual measurements of the  
134 same vial containing the dissolution solution. Between each individual measurement, the  
135 vial was returned to the hot cell line before a new transfer back to the KED in order to take  
136 into account the uncertainties due to the vial position in the KED. The concentration  
137 determined using the KED system was in g L<sup>-1</sup>. In order to obtain the uranium mass fraction

138 in the pellet, the dissolution solution density is necessary and was performed with an  
139 Anton-Paar DMA 58 density meter.

140 Results evaluation and uncertainties estimation

141 *Results evaluation*

142 The bias was calculated using Eq. (2).

$$\text{Bias (\%)} = \frac{[\text{U}] - \text{cert}}{\text{cert}} \cdot 100 \quad (2)$$

143 Where [U] is the uranium concentration determined by ID-TIMS or KED and *cert* is the  
144 ILC assigned value.

145 The precision was evaluated by computing the Relative Standard Deviation (RSD) of all  
146 the measurements.

147 Eq. (3) was used to determine if the analytical method has a statistically significant bias  
148 [15]. If the normalized bias (*NB*) was lower than 2, the method was considered having no  
149 statistically significant bias.

$$\text{NB} = \frac{|[\text{U}] - \text{cert}|}{\sqrt{u^2([\text{U}]) + u_{\text{cert}}^2}} \quad (3)$$

150 Where  $u([\text{U}])$  is the uranium concentration uncertainty determined by ID-TIMS or KED at  
151  $k = 1$  and  $u_{\text{cert}}$  is the assigned value uncertainty at  $k = 1$ .

152

153

## 154 Uranium concentration uncertainty using ID-TIMS

155 The uranium concentration uncertainty ( $u([U]_s)$ , at  $k = 1$ ), determined by ID-TIMS shown  
 156 in Eq. (4), was estimated by combining the uncertainties from each term of Eq. (1). The  
 157 terms of Eq. (1) were considered as not correlated.

$$\begin{aligned}
 u^2([U]_s) = & \left( \frac{m_T}{m_S} \cdot \frac{(235)_T}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{R_{\text{Mix}} - R_T}{1 - R_{\text{Mix}} \cdot R_S} \right)^2 \cdot u^2([U]_T) \\
 & + \left( [U]_T \cdot \frac{1}{m_S} \cdot \frac{(235)_T}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{R_{\text{Mix}} - R_T}{1 - R_{\text{Mix}} \cdot R_S} \right)^2 \cdot u^2(m_T) \\
 & + \left( -[U]_T \cdot \frac{m_T}{m_S^2} \cdot \frac{(235)_T}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{R_{\text{Mix}} - R_T}{1 - R_{\text{Mix}} \cdot R_S} \right)^2 \cdot u^2(m_S) \\
 & + \left( [U]_T \cdot \frac{m_T}{m_S} \cdot \frac{1}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{R_{\text{Mix}} - R_T}{1 - R_{\text{Mix}} \cdot R_S} \right)^2 \cdot u^2((235)_T) \\
 & + \left( -[U]_T \cdot \frac{m_T}{m_S} \cdot \frac{(235)_T}{(238)_S^2} \cdot \frac{M_S}{M_T} \cdot \frac{R_{\text{Mix}} - R_T}{1 - R_{\text{Mix}} \cdot R_S} \right)^2 \cdot u^2((238)_S) \\
 & + \left( [U]_T \cdot \frac{m_T}{m_S} \times \frac{(235)_T}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{(R_{\text{Mix}} - R_T) \cdot R_{\text{Mix}}}{(1 - R_{\text{Mix}} \cdot R_S)^2} \right)^2 \cdot u^2(R_S) \\
 & + \left( -[U]_T \cdot \frac{m_T}{m_S} \cdot \frac{(235)_T}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{1}{1 - R_{\text{Mix}} \cdot R_S} \right)^2 \cdot u^2(R_T) \\
 & + \left( [U]_T \cdot \frac{m_T}{m_S} \cdot \frac{(235)_T}{(238)_S} \cdot \frac{M_S}{M_T} \cdot \frac{1 - R_T \cdot R_S}{(1 - R_S \cdot R_{\text{Mix}})^2} \right)^2 \cdot u^2(R_{\text{Mix}}) + u^2(\text{dilution}) \\
 & + u^2(\text{precision})
 \end{aligned} \tag{4}$$

158 In Eq. (4), the molar masses uncertainties were considered negligible and were not taken  
 159 into account for the final uncertainty computation. The term  $u(\text{precision})$  is the method  
 160 precision, corresponding to the standard deviation of the 3 groups of determination. The  
 161 term  $u(\text{dilution})$  is the sample dilution uncertainty. The concentration ( $[U]_T$ ), isotope ratio  
 162 ( $R_T$ ) and isotope abundance ( $(235)_T$ ) related to the spike solution and their respective  
 163 uncertainties were given in the CRM certificate (Table 1). The sample isotope ratio ( $R_S$ )  
 164 and isotope abundance ( $(238)_S$ ) were estimated in a previous study by the TE method for

165 the  $^{235}\text{U}/^{238}\text{U}$  isotope ratio and the classical method with multi-dynamic sequences for the  
166  $^{234}\text{U}/^{238}\text{U}$  and  $^{236}\text{U}/^{238}\text{U}$  isotope ratios (required to compute the  $^{238}\text{U}$  isotope abundance)  
167 [10].

168 The sample ( $m_S$ ) and spike ( $m_T$ ) masses uncertainties involved in the mixture were  
169 estimated by taking into account the trueness, drift and resolution of the scale and the  
170 method precision (Eq. (5)) [16]. The scale trueness and drift were estimated with a type B  
171 evaluation and a uniform law (first term in Eq. (5)). This term corresponds to the maximum  
172 permissible errors ( $MPE$ ) fixed by the laboratory divided by square root of 3. The scale  
173 resolution ( $q$ ) uncertainty was also estimated using a type B evaluation and a uniform law  
174 (second term in Eq. (5)). The method precision, including scale precision and operator  
175 effect, corresponds to the experimental standard deviation obtained from a series of  $n$   
176 measurements divided by the square root of  $n$  and was a type A evaluation (third term in  
177 Eq. (5)). The air buoyancy (Archimedes' principle) was taken into account for the masses  
178 correction and its uncertainty was considered negligible.

$$u^2(m) = \left(\frac{MPE}{\sqrt{3}}\right)^2 + \left(\sqrt{\frac{q^2}{12}}\right)^2 + \left(\frac{s}{\sqrt{n}}\right)^2 \quad (5)$$

179 The mixture isotope ratio uncertainty was estimated at 0.14 % ( $k = 2$ ), corresponding to  
180 the International Target Value (ITV) for a  $^{235}\text{U}/^{238}\text{U}$  isotope ratio of about 1 [17]. A similar  
181 uncertainty was previously estimated for a  $^{235}\text{U}/^{238}\text{U}$  isotope ratio close to 1 in these  
182 working condition using the TE method with FC 11 [14].

### 183 *Investigation on the optimum isotope ratio of the mixture for the* 184 *ID-TIMS*

185 The ID-TIMS require an optimum mixture isotope ratio ( $R_{Mix}$ ), which can be determined  
186 theoretically [18]. The ID-TIMS parameters other than the mixture isotope ratio ( $[\text{U}]_T$ ,  $m_T$ ,  
187  $m_S$ ,  $(235)_T$ ,  $(238)_S$ ,  $R_S$  or  $R_T$ ) depend only on the sample and spike characteristics and on

188 the mixture masses. These parameters were assumed constant for the theoretical study  
189 hereafter described and equal to the values presented in Table 1.

190 Using Eq. (1) and the fixed parameters given in the Table 1, it is possible to compute using  
191 a spreadsheet the sample concentration for a given mixture isotope ratio. This  $^{238}\text{U}/^{235}\text{U}$   
192 mixture isotope ratio is between the sample (0.0581) and spike (50.87) isotope ratios.

193 Likewise, the sample concentration uncertainty can be computed using Eq. (4) for each  
194 mixture isotope ratio included between the  $^{238}\text{U}/^{235}\text{U}$  sample and spike isotope ratios. The  
195  $u(\textit{precision})$  and  $u(\textit{dilution})$  terms were assumed to be constant in this theoretical study.  
196 The  $u(\textit{precision})$  term was set at 0.06 % ( $k = 1$ ), which is a value previously obtained in  
197 our laboratory for actinide concentration determination by ID-TIMS [14]. The  $u(\textit{dilution})$   
198 term was set at 0.034 %, which is an arbitrary value commonly observed in our laboratory  
199 for this type of dilution. These 2 set values are not function of the mixture isotope ratio and  
200 do not change the optimum determination. By combining Eq. (1) and (4), the concentration  
201 relative uncertainty ( $u([\text{U}]/[\text{U}], k = 1)$ ), as a function of the mixture isotope ratio, can be  
202 drawn from the spreadsheet (Fig. 1). The optimum mixture isotope ratio is obtained when  
203 the function reaches its minimum.

204 In our analysis conditions, the optimum mixture isotope ratio is obtained for a  $^{238}\text{U}/^{235}\text{U}$   
205 isotope ratio of 1.7 (Fig. 1). Experimentally, it is possible to adjust the sample and spike  
206 solution amount in the mixture to be closer to the optimum mixture isotope ratio. In  
207 addition, this theoretical study shows that the sample concentration relative uncertainty is  
208 almost stable around the optimum isotope ratio. It shows that the uncertainty increases  
209 large when the mixture isotope ratio approaches the spike or sample isotope ratios. The  
210 mixture isotope ratios from 1.08 to 1.23 (which are the ratios measured in this study),  
211 increase the concentration uncertainty only by 0.09 % to 0.18 % compared to the best  
212 theoretically achievable uncertainty.

213 *Uranium concentration uncertainty using KED analysis*

214 The main parameters to take into account for the concentration uncertainty by KED  
215 measurement are: the instrument variability (X-ray tube, detector), the vial position, the  
216 counting time precision, the vial length, the sample matrix and the theoretical values of  
217 mass attenuation coefficients. The uranium concentration relative uncertainty  $u([U])/[U]$   
218 ( $k = 1$ ) was estimated using Eq. (6). The first three terms correspond to the KED  
219 measurement uncertainty. The fourth term is the density uncertainty of the dissolution  
220 solution at  $k = 1$  and also includes the density difference between the density and KED  
221 measurements as the two measurements were not performed at the exact same temperature.  
222 The first term of Eq. (6) takes into account the measurement precision component. It  
223 includes the uncertainty from random effects (instrument variability, vial position and  
224 counting time precision) and is given by the RSD of all the measurements. The second and  
225 third terms take into account the systematic error components (vial length, sample matrix  
226 and theoretical values of mass attenuation coefficients). The measurement trueness  
227 ( $u^2(\text{trueness})/(\text{trueness})^2$ ) is calculated with Eq. (7). The CRM were chosen so that they  
228 respected different criteria: same element, similar analyte concentration and matrix. The  
229 EQRAIN solutions, provided by the CEA/CETAMA during ILC were used [4]. They were  
230 chosen because: the analyte is the same, the uranium concentrations are close (about 250  
231 and 360 g L<sup>-1</sup> for the EQRAIN solutions and about 280 g L<sup>-1</sup> for the dissolution sample)  
232 and the matrix is the same (uranium only in HNO<sub>3</sub> media). Only the EQRAIN solutions  
233 analyzed using the same configuration and the same detector as the dissolution sample (as  
234 the detector was changed in 2016) were used to compute the trueness term.

235 The term  $u^2(\text{ref})/(\text{ref})^2$  corresponds to the CRM (EQRAIN U) concentration relative  
236 uncertainty at  $k = 1$ .

$$\frac{u^2([U])}{([U])^2} = \frac{u^2(\bar{x})}{\bar{x}^2} + \frac{u^2(\text{trueness})}{(\text{trueness})^2} + \frac{u^2(\text{ref})}{(\text{ref})^2} + \frac{u^2(\text{density})}{(\text{density})^2} \quad (6)$$

$$\frac{u(\text{trueness})}{\text{trueness}} = \frac{\text{Maximum bias on CRM}}{\sqrt{3}} \quad (7)$$

237        **Results and discussion**

238        ID-TIMS

239        The results are presented in Table 2. The RSD of the 3 groups of determinations is equal  
240        to 0.12 %. The result obtained for the uranium mass fraction of the pellet is  
241         $0.8835 \pm 0.0028 \text{ g g}^{-1}$  ( $k = 2$ ). The bias is equal to 0.30 %. The relative uncertainty,  
242        computed with Eq. (4), is estimated at 0.32 % and is in compliance with the safeguard ITV  
243        (0.36 %,  $k = 2$ ) [17]. The relative contribution study (expressed in percent) of the main  
244        uncertainty sources for the uranium mass fraction is given in Fig. 2. The spike  
245        concentration uncertainty (0.04 %) has a very low contribution and is a minor source of  
246        uncertainty (4.2 % of the total budget). The main uncertainty sources are the precision  
247        (59.4 % of the total budget), the dilution to obtain a sample concentration of about  $10 \mu\text{g g}^{-1}$   
248        ( $^{238}\text{U}/^{235}\text{U}$  sample - spike mixture isotope ratio  
249        determination (14.1 % of the total budget).

250        KED

251        A uranium concentration of  $278.4 \pm 2.1 \text{ g L}^{-1}$  ( $k = 2$ ) was found for the dissolution solution.  
252        Considering the dissolution solution density and the mass of the initial pellet before  
253        dissolution, the uranium mass fraction of the pellet is estimated to  $0.8807 \pm 0.0067 \text{ g g}^{-1}$   
254        ( $k = 2$ , Table 2). The RSD of the 5 individual determinations is 0.33 %. The bias is below  
255        0.1 %. The uranium mass fraction determined using the KED system is in compliance with  
256        the assigned value provided by the IAEA: the normalized bias is equal to 0.09. The relative  
257        uncertainty is estimated to 0.76 % and is slightly higher than the ITV (0.56 %,  $k = 2$ ) [17].  
258        The study of the main uncertainty sources shows the total uncertainty is mainly due to the  
259        precision (76 % of the total uncertainty). The trueness component (18 % of the total  
260        uncertainty) and the CRM reference uncertainty component (6 % of the total uncertainty)  
261        have a minor contribution to the total uncertainty. The density uncertainty contribution is  
262        negligible: below 0.2 % of the total uncertainty.

263 ID-TIMS and KED comparison

264 The ID-TIMS and the KED perform uranium assay with high accuracy. However, in terms  
265 of analytical performance, implementation, analysis time and generated radioactive  
266 effluent, they have different assets (Table 3).

267 The KED requires high uranium concentration solution in order to obtain the lowest  
268 uncertainties: uranium concentration higher than  $50 \text{ g L}^{-1}$  is needed to get relative expanded  
269 uncertainties ( $k = 2$ ) lower than 1%. The TIMS analyses require a lower uranium amount  
270 than the KED: less than  $1 \mu\text{g}$  of uranium is required for TIMS analyses. Solutions  
271 containing less than  $1 \text{ mg L}^{-1}$  of uranium can be analyzed by ID-TIMS and the ITV can  
272 still be achieved (0.36 %,  $k = 2$ ). The ID-TIMS gives directly the uranium mass fraction in  
273 the nuclear materials, whereas the density determination for the KED system is required  
274 beforehand. Both techniques showed no significant bias compared to the assigned value:  
275 the normalized biases are below 2. The ID-TIMS uncertainty estimation (0.32 %,  $k = 2$ )  
276 was in compliance with the ITV and was lower than the KED uncertainty estimation  
277 (0.76 %,  $k = 2$ ).

278 Both techniques require an identical pellet dissolution step. KED analysis can be performed  
279 directly on the dissolution solution (matrix  $8 \text{ mol L}^{-1} \text{ HNO}_3$ ). No solution treatment, no  
280 dilution and no calibration are required [4]. Cross contamination is then limited. As the  
281 KED system is a nondestructive technique, the same aliquot can be used to measure  
282 uranium with the KED system and be used for others analyses like isotopic, radiometric or  
283 density measurements (in some situations, taking into account evaporation loss is  
284 necessary). High activity solution can be analyzed with the KED system as it is installed  
285 in a shielded line. The TIMS can be installed in glove box environment but not in a hot  
286 cell: high activity solutions require a dilution step in order to be suitable for working in  
287 glove box. The dilution step is also required in order not to consume too much of the  
288 solution spike. Indeed, for spiking the dissolution solution directly would have a  
289 prohibitive cost as the spike solution quantity would be important. For very irradiating  
290 samples such as high active liquid waste, the spiking step must be done in the hot cell  
291 before the sample transfers.

292 As the analytical protocol is very different, both techniques present different analysis time.  
293 For KED, about 5 minutes are required to prepare the sample in the analytical vial. The  
294 analysis time is about 1 h 15 for 5 analysis of the same sample (about 15 min for 1  
295 determination). Considering the fast analysis time of the KED system; the dissolution  
296 solution density determination is the most critical step to the total analysis time (about  
297 5 h 20). The ID-TIMS analysis overall process takes longer than the KED analysis. The  
298 dilution step is performed in 1 h. The sample preparation, including sample - spike mixture  
299 preparation, sample - spike mixture evaporation, deposit onto the filament, filament  
300 loading inside the TIMS source take about 8 h in total. Then, the vacuum must be drawn  
301 inside the source (about 4 h). This step is generally not critical in terms of operator working  
302 time as the filament loading can be done at the end of a working day and the vacuum inside  
303 the source can be obtain overnight. The analysis time of the 10 sample - spike mixtures  
304 with the TE method requires about 9 h of work: Faraday cups calibration duration is about  
305 10 min per day of use, heating filament about 10 min per analysis, signal optimization  
306 about 5 min per analysis, the analysis time duration is about 30 min and the shutdown is  
307 about 2 min per analysis. The ID-TIMS methodology requires also the sample isotope  
308 determination. In total, we have calculated the analysis time for the sample isotope  
309 determination with TIMS at 19 h [10]. The total ID-TIMS analysis time is estimated at  
310 41 h, thus about 8 times higher for the ID-TIMS than KED.

311 The KED system does not consume the sample solution as it is a nondestructive technique  
312 and generates no effluent as there is no dilution step. However, this technique consumes  
313 sample solution indirectly as the density determination requires about 2 mL that cannot be  
314 reused due to possible cross-contamination. The effluent generated for cleaning the density  
315 meter represents about 10 mL. The total effluent generated was 12 mL, only due to density  
316 determination. Moreover, the used KED system requires at least 1.5 mL of solution, which  
317 can be a disadvantage when the available sample volume is low. The TIMS analyses need  
318 a dilution step as uranium concentration is high (higher than  $1 \text{ g.L}^{-1}$ ). In this study, 600  $\mu\text{L}$   
319 of the dissolution solution was used, corresponding to  $3 \times 200 \mu\text{L}$  for the 3 determinations  
320 by ID-TIMS. 200  $\mu\text{L}$  of solution (about 300 mg of uranium) is enough to obtain an accurate  
321 weight without generating too much nuclear effluent. The generated effluent was estimated  
322 at 60 mL for one determination. In this study, the 3 determinations generated 180 mL of

323 effluent, only due to the dilution step. The ID-TIMS (0.6 mL) consumes less sample than  
324 the KED system (2 mL required only for the density determination). However, the KED  
325 system generates 15 times less volume of radioactive effluent.

## 326 **Conclusions**

327 This study presents a comparison between the ID-TIMS and the KED system for uranium  
328 concentration determination with high accuracy. The results obtained with both techniques  
329 are in compliance with the assigned value. Both techniques can be used for metrological  
330 analysis of high uranium concentration while having different assets. The minimal  
331 measurable concentration with ID-TIMS is lower than the KED's. The analysis time is  
332 about 8 time higher for the ID-TIMS than for the KED. The ID-TIMS generates about 15  
333 times more radioactive effluent. However, in our studies, only the estimated uncertainties  
334 obtained with ID-TIMS are in compliance with the ITV. The estimated uncertainties  
335 obtained with the KED are slightly higher.

336 The main approach to decrease the KED uncertainty, in order to be closer to the ITV, is to  
337 improve the measurement precision as it is the main uncertainty source (76 % of the total  
338 uncertainty budget). The vial position which may vary slightly from an analysis to another  
339 in the KED system and the instrument variability (X-ray tube, detector) are the main  
340 parameters influencing the method precision as well as the final uncertainty. An approach  
341 to decrease the precision due to the instrument variability would be increasing the  
342 measurement time. Even with a measurement time of 1 hour (15 min at present), the KED  
343 would be faster than ID-TIMS. The precision due to the instrument variability would be  
344 improved. However the precision due to the vial position would be the same whatever the  
345 measurement time. The method used in the laboratory is a calibration free method, which  
346 slightly increases the estimated uncertainties compared to methods using standard  
347 calibration. A calibration would correct all physical phenomena including the ones that are  
348 difficult to take into account in theoretical models. However, this correction would only  
349 modify the method trueness. Therefore this would not change much the final uncertainty  
350 as the trueness has a low contribution part in the final uncertainty (18 % of the total

351 uncertainty). The KED data were also processed using the algorithms with fitting windows  
352 as proposed in [5] in order to show the influence on the repeatability. The repeatability  
353 were similar for either the algorithms using fitting windows (RSD = 0.47 %) or our in-  
354 house algorithm (RSD = 0.33 %). In addition, the KED system used for this study is  
355 designed for dissolution or separation R&D studies that need quick results during the  
356 experiments. Our KED system, implemented in hot-cell facility and used as a calibration  
357 free method, is fully compatible for this type of application.

### 358 **Acknowledgements**

359 We are grateful to Dr. S. Baghdadi (IRSN/PSN-EXP/SSRD/BTE) for her precious advice  
360 on the present paper. We would like to thank all of the IAEA staff (Seibersdorf-Austria)  
361 who organized the interlaboratory comparison and all the KED staff for the KED  
362 measurements.

363 **Table**

364 Table 1: Parameters used for the theoretical study of our studied case in order to find the  
 365 optimum sample – spike mixture isotope ratio in the ID-TIMS

Parameters	Value	Uncertainty (k=1)
$[U]_T$	9.5044	0.0039
$m_T$	0.246555	0.000059
$m_S$	0.218465	0.000059
$M_S$	237.992	-
$M_T$	235.201	-
$(^{235})_T$	93.1760	0.0034
$(^{238})_S$	98.053	0.010
$R_S$	0.019658	0.000015
$R_T$	0.058065	0.000017
Parameters	Value	Relative Uncertainty (k=1)
$R_{Mix}$	Between 0.0586 and 50.87	0.07 %

366

367

368

369

370

371

372 Table 2: Uranium mass fraction determined by ID-TIMS and KED obtained on the  
 373 uranium pellet.  $U_{ass\_val}$  corresponds to the assigned value uncertainty ( $k = 2$ ) and  $U$   
 374 corresponds to the total uncertainty ( $k = 2$ ). NB corresponds to the normalized bias

Method	Uranium mass fraction	
	ID-TIMS	KED
Assigned value ( $\text{g g}^{-1}$ )	0.88110	0.88110
Results ( $\text{g g}^{-1}$ )	0.8835	0.8807
<i>RSD</i> (%)	0.12	0.33
<i>NB</i>	1.8	0.09
Bias (%)	0.28	-0.04
$U_{ass\_val}$ (% , $k = 2$ )	0.05	0.05
$U$ (% , $k = 2$ )	0.32	0.76
<i>ITV</i> (% , $k = 2$ )	0.36	0.56


375

376 Table 3: Experimental duration and effluents

Technique	ID-TIMS	KED
		<u>Duration</u>
Dilution step	1 h	/
Sample preparation	12 h	5 min
Analysis	9 h	1 h 15 ( $5 \times 15$ min)
Density	/	4h
Isotope determination	19 h	/
Total	41 h	5h20
		<u>Waste management</u>
Sample solution consumption	600 $\mu\text{L}$	- 2 mL for the density determination - 0 mL for the KED analysis
Generated effluent	180 mL for the dilution step	10 mL for the density determination

377


378 **Figure**

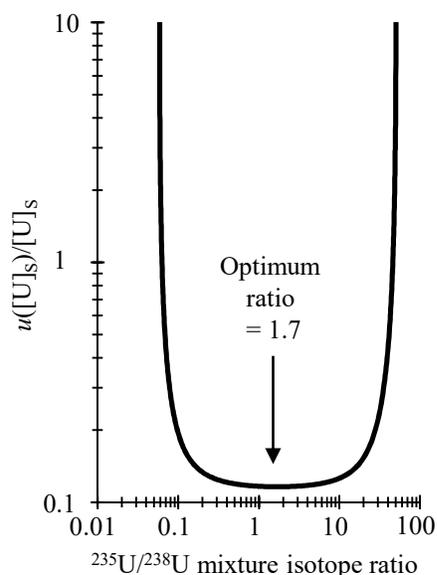

379

380 **Fig. 1:** Schematic of the uranium concentration relative uncertainty ( $u([U]_s)/[U]_s$ , %,  $k = 1$ ) as a function of the mixture isotope ratio in order to find the optimum mixture  
 381 isotope ratio in our studied case. The  $[U]_s$  and  $u([U]_s)$  terms were calculated using Eq.  
 382 (1) and Eq. (4), respectively, considering all the isotope dilution terms constants (Table 1)  
 383 except the mixture isotope ratio value  
 384

385

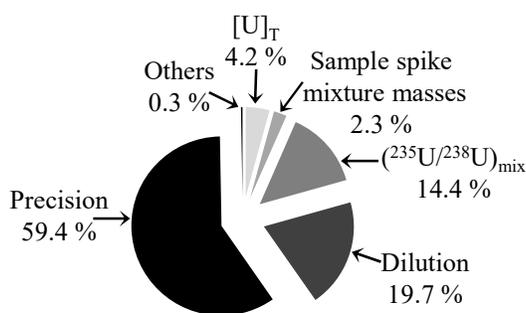

386

387 **Fig. 2:** Relative contribution (%) of the main uncertainty sources for the uranium mass  
 388 fraction determined by ID-TIMS

389

390        **References**

- 391    1.    Godoy MLDP, Godoy JM, Roldão LA, Tauhata L (2009) Determination of total  
392        content and isotopic compositions of plutonium and uranium in environmental  
393        samples for safeguards purposes by ICP-QMS. *J Environ Radioact* 100:613–625.  
394        doi: 10.1016/j.jenvrad.2009.04.012
- 395    2.    Metzger SC, Ticknor BW, Rogers KT, Bostick DA, McBay EH, Hexel CR (2018)  
396        Automated Separation of Uranium and Plutonium from Environmental Swipe  
397        Samples for Multiple Collector Inductively Coupled Plasma Mass Spectrometry.  
398        *Anal Chem* 90:9441–9448. doi: 10.1021/acs.analchem.8b02095
- 399    3.    Boulyga S, Konegger-Kappel S, Richter S, Sangély L (2015) Mass spectrometric  
400        analysis for nuclear safeguards. *J Anal At Spectrom* 30:1469–1489. doi:  
401        10.1039/c4ja00491d
- 402    4.    Esbelin E, Rivier C, Duhamel G (2018) High concentration measurements of U and  
403        Pu with nondestructive and standardless K-edge densitometer device. *X-Ray*  
404        *Spectrom* 47:22–33. doi: 10.1002/xrs.2802
- 405    5.    Ottmar H, Eberle H (1991) The hybrid K-edge/K-XRF densitometer: Principles-  
406        design-performance. KFK Rep. 4590, Karlsruhe
- 407    6.    Quemet A, Maillard C, Ruas A (2015) Determination of zirconium isotope  
408        composition and concentration for nuclear sample analysis using Thermal Ionization  
409        Mass Spectrometry. *Int J Mass Spectrom* 392:34–40. doi:  
410        10.1016/j.ijms.2015.08.023
- 411    7.    Aggarwal SK (2016) Thermal ionisation mass spectrometry (TIMS) in nuclear  
412        science and technology – a review. *Anal Methods* 8:942–957. doi:  
413        10.1039/c5ay02816g
- 414    8.    Quemet A, Maloubier M, Dalier V, Ruas A (2014) Development of an analysis  
415        method of minor uranium isotope ratio measurements using electron multipliers in

- 416 thermal ionization mass spectrometry. *Int J Mass Spectrom* 374:26–32. doi:  
417 10.1016/j.ijms.2014.10.008
- 418 9. Quemet A, Maloubier M, Ruas A (2016) Contribution of the Faraday cup coupled  
419 to 1012 ohms current amplifier to uranium 235/238 and 234/238 isotope ratio  
420 measurements by Thermal Ionization Mass Spectrometry. *Int J Mass Spectrom*  
421 404:35–39. doi: 10.1016/j.ijms.2016.04.005
- 422 10. Quemet A, Ruas A, Dalier V, Rivier C (2019) Development and comparison of high  
423 accuracy thermal ionization methods for uranium isotope ratios determination in  
424 nuclear fuel. *Int J Mass Spectrom* 438:166–174. doi: 10.1016/j.ijms.2019.01.008
- 425 11. Richter S, Kühn H, Aregbe Y, Hedberg M, Horta-Domenech J, Mayer K, Zuleger  
426 E, Bürger S, Boulyga S, Köpf A, Poths J, Mathew K (2011) Improvements in routine  
427 uranium isotope ratio measurements using the modified total evaporation method  
428 for multi-collector thermal ionization mass spectrometry. *J Anal At Spectrom*  
429 26:550–564. doi: 10.1039/C0JA00173B
- 430 12. Wegener MR, Mathew KJ, Hasozbek A (2013) The direct total evaporation (DTE)  
431 method for TIMS analysis. *J Radioanal Nucl Chem* 296:441–445. doi:  
432 10.1007/s10967-012-2182-3
- 433 13. Mathew KJ, O'Connor G, Hasozbek A, Kraiem M (2013) Total evaporation method  
434 for uranium isotope-amount ratio measurements. *J Anal At Spectrom* 28:866–876.
- 435 14. Quemet A, Ruas A, Dalier V, Rivier C (2018) Americium isotope analysis by  
436 Thermal Ionization Mass Spectrometry using the Total Evaporation Method. *Int J*  
437 *Mass Spectrom* 431:8–14. doi: 10.1016/j.ijms.2018.05.017
- 438 15. AFNOR (2009) NF T90-210 norm : Water quality - Protocol for the initial method  
439 performance assesment in a laboratory.
- 440 16. Désenfant M, Priel M, Rivier C (2005) Evaluation des incertitudes des résultats  
441 d'analyse. Ref P105 V1. Les Tech l'Ingénieur 1–17.

- 442 17. International Atomic Energy Agency (2010) International Target Values 2010 for  
443 Measurement Uncertainties in Safeguarding Nuclear Materials - STR368. Vienna,  
444 Austria
- 445 18. Irrgeher J, Vogl J, Santner J, Prohaska T (2015) Chapter 8 Measurement Strategies.  
446 In: Sect. F. Mass Spectrom. Elem. Isot. Anal. The Royal Society of Chemistry, pp  
447 126–151
- 448